

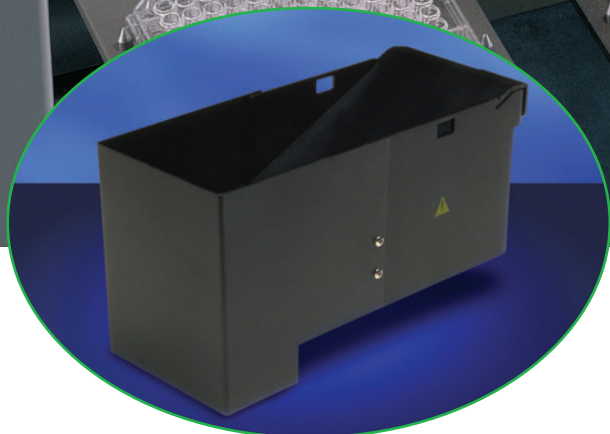


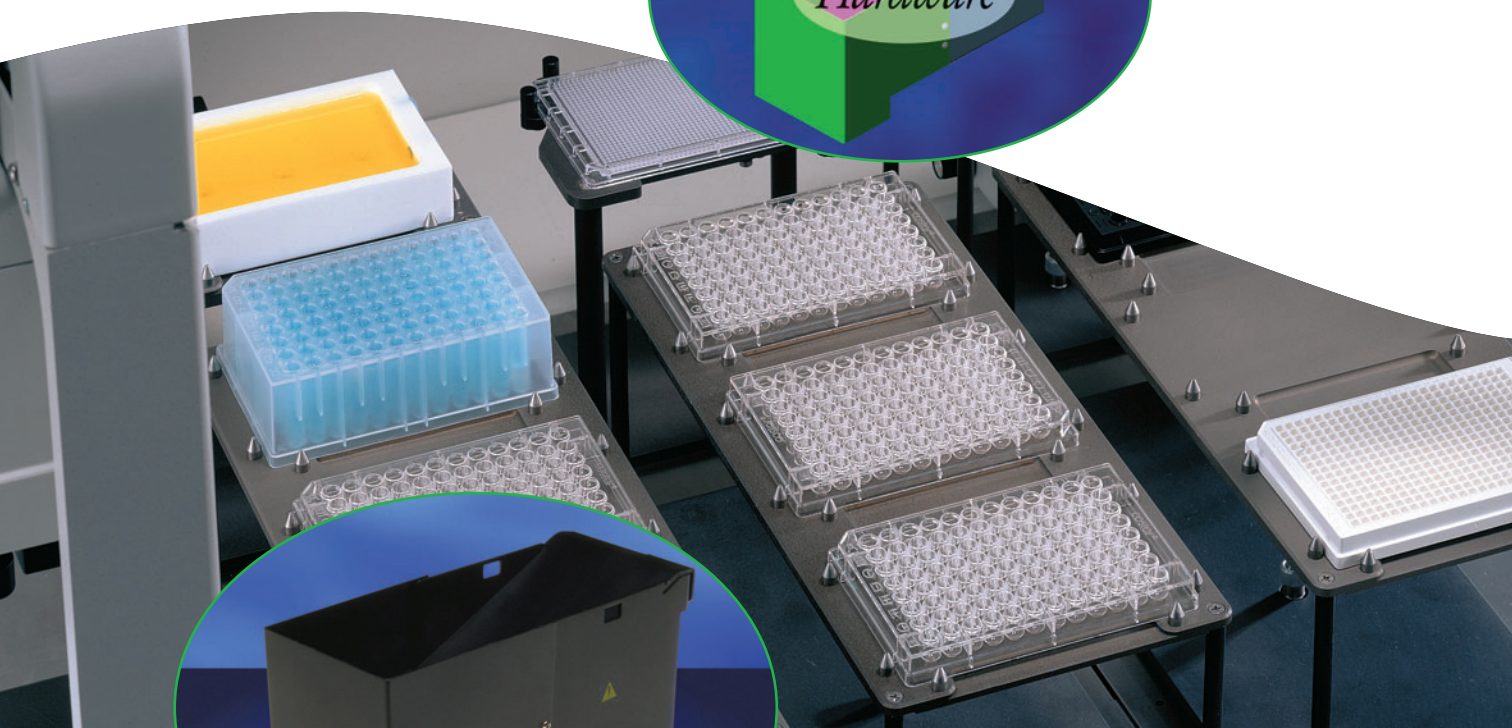


## Automated Labware Positioning.

### BIOMEK® FX

STATIC ALPs AND  
WASTE ALP

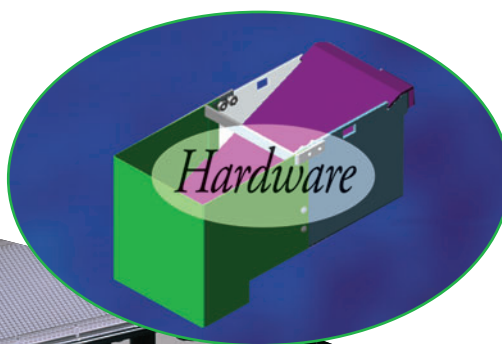

*Biomek FX Waste ALP*

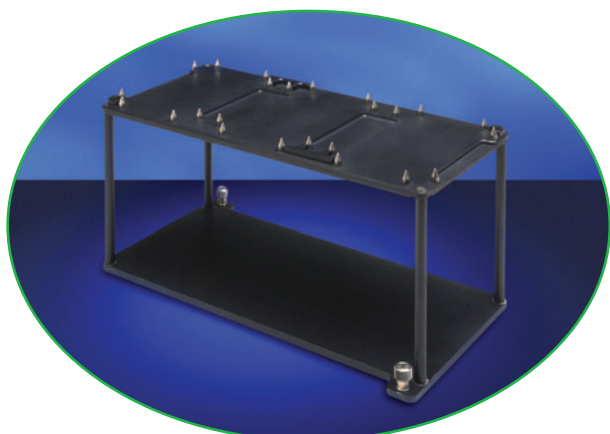

*Biomek FX 3-position Static ALP*


The Biomek FX has a large flexible worksurface. Integral to the system are Automated Labware Positioners (ALPs) that are uniquely designed to fit seamlessly onto the deck. These ALPs perform application-specific tasks and make it simple to configure the system for a wide variety of applications.

- Waste ALP can be positioned on the outer columns of the deck for disposal of expired labware, tips and tip boxes.
- A variety of Static ALPs can be used in combination with active ALPs to create application-specific decks.


ALPs, one more reason the Biomek FX is the *smart solution* for your liquid handling needs.

## Automated Labware Positioning.

Automated Labware Positioners (ALPs) are user interchangeable for easy reconfiguration, enhanced productivity and reliability.

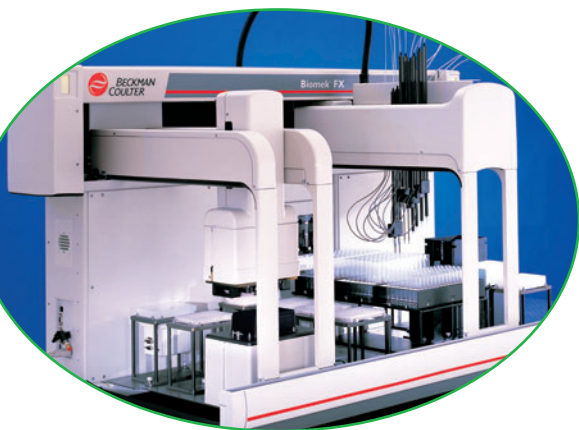

Biomek FX 1x3 ALP


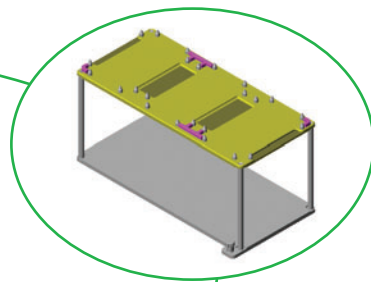


### Static ALPs

- Useful for locating a wide variety of labware that conforms to SBS footprint standards
- Choices include 1x1, 1x3, 4x3 and 4x4 for labware that will be moved by the gripper
- A 1x5 is available for labware such as reagent reservoirs that will not be moved
- Combine with active ALPS to create application-specific deck layout

### Waste ALP

- Useful for automated disposal of expired labware, tips and tip boxes
- Use multiple Waste ALPs to segregate waste streams
- Waste ALP required on all Span-8 systems for disposable tip removal

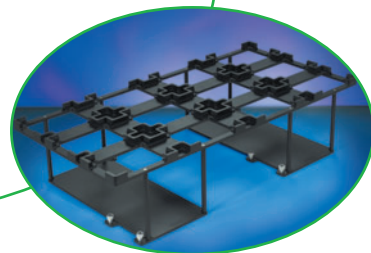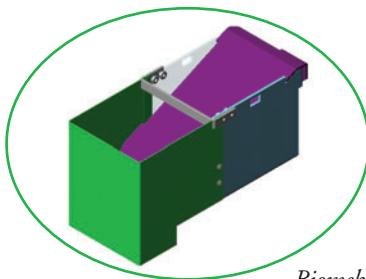

Biomek FX 4x3 ALP

Biomek FX Waste ALP

## Ordering Information

719357 Standard Single-Position ALP

719358 Standard Three-Position ALP

719948 4x3 Static ALP Kit

719360 High-Density Sixteen-Position ALP

*The above four ALPs are for positioning standard 96 and/or 384 well microtiter plates as singles or in stacks with spacing to support Biomek FX gripper and/or robotic manipulation.*

719359 Standard Five-Position ALP

*A five-position ALP for positioning of five (5) reservoirs and/or microplates as singles. The five-position ALP DOES NOT support Biomek FX gripper and/or robotic manipulation.*

719347 Waste ALP (Multichannel Pod)

719590 Waste ALP (Span-8 Pod)

*Includes waste ALP, chute, and disposable bags. The waste ALP chute is for the disposal of tips, full or empty tip racks and microtiter plates from the work area of the Biomek FX.*



[www.beckmancoulter.com/smartsolutions](http://www.beckmancoulter.com/smartsolutions)

*Developing innovative solutions in genetic analysis, drug discovery, and instrument systems.*

**Innovate** Automate  
SIMPLIFY



Beckman Coulter, Inc. • 4300 N. Harbor Boulevard, Box 3100 • Fullerton, California 92834-3100

Sales: 1-800-742-2345 • Service: 1-800-551-1150 • Telex: 678413 • Fax: 1-800-643-4366 • [www.beckmancoulter.com](http://www.beckmancoulter.com)

### Worldwide Life Science Research Division Offices:

**Australia** (61) 2 9844-6000 **Canada** (905) 819-1234 **China** (86) 10 6515 6028 **Eastern Europe, Middle East, Africa** (41) 22 994 07 07

**France** 01 49 90 90 00 **Germany** (89) 35870-0 **Hong Kong** (852) 2814 7431/2814 0481 **Italy** 02-953921 **Japan** 03-5404-8359

**Mexico** 525-605-77-70 **Netherlands** 0297-230630 **Singapore** (65) 6339 3633 **South Africa, Sub-Saharan Africa** (27) 11-805-2014/5 **Spain** (34) 91 3836080

**Sweden** 08-564 85 900 **Switzerland** 0800 850 810 **Taiwan** (886) 2 2378 3456 **Turkey** 90 216 309 1900 **U.K.** 01494 441181 **U.S.A.** 1-800-742-2345

B2002-5057-CB-10

© 2002 Beckman Coulter, Inc.

Printed in U.S.A. on recycled paper.